



Instructions for Applying the **NANRIC** HE0 - 2 Heel Extensions

NANRIC *State of the Art Treatment for Foals with Weak Flexor Tendons*




FOOT PREPARATION

- Trim the foot removing the center of pivot. This will offer support to weak flexor tendons.
- Round the ground surface of the hoof to accommodate the seating of the foot into the cuff.
- Roughen the wall slightly with the smooth side of the rasp to provide an adequate glue contact surface.
- Thoroughly clean the dirt from all the hoof cracks and cover potential abscesses with a small piece of tape or Keratex Putty to prevent sealing any bacteria under the glue.
- Fit the shoe to the foot, and trim any area of the cuff that touches the coronary band.
- Use a marker to make a matching mark on the foot and shoe so you can quickly line up the marks when gluing.

SHOE APPLICATION

- If you are working with an animal that is fractious, you may have your vet administer a sedative.
- Clean the hoof wall with denatured alcohol and ensure the foot is dry. Use a heat gun if necessary.
- Open the **NANRIC** Mega Bond™ package and unscrew the locking collar and remove the storage cap. Attach a mixing tip by pressing on to the open cartridge. Secure the locking collar over the top of the attached tip with a quarter turn. Insert the cartridge into the dispensing gun and dispense a small amount through the mixing tip to equalize the cartridge.
- Scan the QR code on the back of the cartridge for more detailed adhesive instructions.
- Apply the shoe to the clean foot.
- Insert the glue tip into the glue ports to apply **NANRIC** Mega Bond™ to the inside of the cuff wall only.
- **Mega Bond™ sets in 30 seconds, so remember to work quickly.** Total cure time is 2-3 minutes. If you need a longer set time, Mega Bond™ can be stored in a refrigerator prior to use.
- To apply a second shoe, simply apply a new tip to the Mega Bond™ cartridge.
- When finished, apply the storage plug to the Mega Bond™. The remaining adhesive can be used at another time.

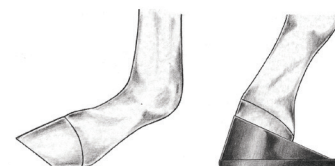

SHOE REMOVAL

- Use nippers or a rasp to separate the entire cuff from the horn wall and then pull the shoe.
- if a reset is necessary, we recommend applying a new heel extension as it is difficult to clean out old glue residue.

RECOMMENDED SHOEING SCHEDULE

- For the best results, follow this shoeing schedule as the foal's feet continue to grow. Exceeding the recommended shoeing time can alter natural foot development.

Age	Maximum Time for Shoes To Be Worn	Days Following Shoe Removal To Re-evaluate
Birth - 2 Weeks	4 Days	7 Days
2-4 Weeks	7-10 Days	7-10 Days
4-8 Weeks	10-14 Days	7-10 Days
8-12 Weeks	14 Days	7-10 Days

SORENESS AFTER APPLICATION

- Fitting the shoe too tightly will pinch the quarters causing slight soreness. If this occurs, pull the shoe and reapply after the foal is sound.
- Sealing bacteria in broken areas along the white line will cause acute lameness. If this occurs, pull the shoe immediately and notify your veterinarian. The abscess will need to be opened and appropriate therapy must be administered.